

10340 Camino Santa Fe, Suite C, San Diego, CA 92121

Tel: (858)412-5988 Fax: (858)368-8716

Technical Supports: <a href="mailto:supports@ixcellsbiotech.com">supports@ixcellsbiotech.com</a>

Orders: orders@ixcellsbiotech.com

## **Product Information**


#### **Human Hepatic Stellate Cells (HHSC)**

Catalog Number	10HU-210	Cell Number	1.0 million cells/vial
Species	Human	Storage Temperature	Liquid nitrogen

# **Description**

Hepatic stellate cells (HSCs) are liver-specific mesenchymal cells, and account for 5~8% of the cells in the liver. HSCs play vital roles in the homeostasis of liver extracellular matrix, repair, regeneration and fibrosis, and control retinol metabolism, storage and release. The stellate cell is the major cell type involved in liver fibrosis in response to liver injury. In healthy liver, HSCs are in a quiescent state, and contains numerous vitamin A lipid droplets, constituting the largest reservoir of vitamin A in the body. When the liver is damaged, HSCs can change into an activated state, which is characterized by proliferation, contractility and chemotaxis. The amount of vitamin A decreases progressively in injured liver. The activated HSCs also secrete collagen scar tissue, which can lead to cirrhosis. In chronic liver disease, prolonged and repeated activation of stellate cells causes liver fibrosis [1,2]. Primary culture of HSC is a valuable tool to study liver fibrosis.

iXCells Biotechnologies provides high quality Human Hepatic Stellate Cells (HHSC), which are isolated from adult human liver and cryopreserved with >1 million cells in each vial. Each vial is sufficient to seed on four 10 cm dishes, and cells will be ready in 5~6 days to be passaged at 1:3 ratio. HUSMC express Desmin and GFAP, and are negative for HIV-1, HBV, HCV, mycoplasma, bacteria, yeast, and fung. HUSMC can further expand for 2-3 population doublings in Stellate Cell Growth Medium (Cat # MD-0014) under the condition suggested by iXCells Biotechnologies.

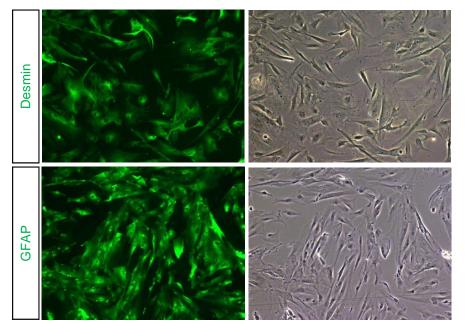

**Figure 1.** Human hepatic stellate cells (HHSC) were stained with HSC markers, including Desmin and GFPA as indicated (green). Phase contrast images of each staining is shown on the right panel.

### **Product Details**

Tissue	Human liver	
Package Size	1.0 million cells/vial	
Shipped	Cryopreserved	
Storage	Liquid nitrogen	
<b>Growth Properties</b>	Adherent	
Media	Stellate Cell Growth Medium (Cat# MD-0014)	

### **Protocols**

#### **Thawing of Frozen Cells**

- 1. Upon receipt of the frozen Human Hepatic Stellate Cells (HHSC), it is recommended to thaw the cells and initiate the culture immediately in order to retain the highest cell viability.
- 2. To thaw the cells, put the vial in 37°C water bath with gentle agitation for ~1 minute until just the tiniest ice crystal is left in the cryovial. Keep the cap out of water to minimize the risk of contamination.
- 3. It is recommended to seed one vial of HHSCs (1.0-1.2 million cells) into four 100mm cell culture dishes in Stellate Cell Growth Medium (Cat # MD-0014).

Note: Do not use cell culture flasks, because HHSCs will become activated faster than culturing in 100mm dishes.

4. Change the medium every other day.

Safety Precaution: it is highly recommended that protective gloves and clothing should be used when handling frozen vials.

#### **Sub-culture Procedure**

- 1. Human Hepatic Stellate Cells (HHSC) can be cultured in Stellate Cell Growth Medium (Cat # MD-0014).
- 2. When cells reach ~80-90% confluence, remove the medium, and wash once with sterile PBS.
- 3. Add ~2.5ml of 0.25% Trypsin-EDTA to the flask and incubate for ~3 minutes at 37°C. Shake and check every minute until cells detach from the bottom of the plate. Neutralize the enzyme by adding 2-3 volumes of cell culture medium.
- 4. Centrifuge 1000rpm (~220g) for 5min and resuspend the cells in desired volume of medium.
- 5. Seed the cells into new culture vessels with 1:4 dilution.
- 6. Change the medium every other day.

#### **TGFb1 Treatment**

- 1. Seed HHSCs into desired culture vessels with Stellate Cell Growth Medium (Cat # MD-0014).
- 2. After HHSCs attached (usually in 24 hours), serum starve cells for 12-24 hours using DMEM basal medium without FBS.

- 3. Treat cells with 2ng TGFb1 in DMEM basal medium for 12-24 hours.
- 4. Wash twice with PBS. The cells will then be ready for downstream analysis.

### References

- [1] Chunyue Yin, et al, J Clin Invest. 2013;123(5):1902–1910. Hepatic stellate cells in liver development, regeneration, and cancer.
- [2] Rockey D. C. . Semin Liver Dis 21(3):337-49. (2001) Hepatic blood flow regulation by stellate cells in normal and injured liver.

#### **Disclaimers**

This product is intended for laboratory research purposes only. It is not intended for use in humans. While iXCells Biotechnologies uses reasonable efforts to include accurate and up-to-date information on this product sheet, we make no warranties or representations as to its accuracy. Citations from scientific literature and patents are provided for informational purposes only. iXCells Biotechnologies does not warrant that such information has been confirmed to be accurate.

This product is sent with the condition that you are responsible for its safe storage, handling, and use. iXCells Biotechnologies is not liable for any damages or injuries arising from receipt and/or use of this product. While reasonable effort is made to insure authenticity and reliability of strains on deposit, iXCells Biotechnologies is not liable for damages arising from the misidentification or misrepresentation of cultures.

© iXCells Biotechnologies 2015. All rights reserved.